



RIGGTEK
Dissolution Test Systems



easyDISS TX8

features and advantages (1)



The easyDISS TX8 is all you really need at a dissolution tester

- Rugged and compact construction for maximum reliability
- Certificate of Compliance for USP and Ph.Eur.
- All relevant parts as standard with serial number and Certificate of Conformance – optional with Certificate of Analysis
- Easy to use and clean
- standard with 8 positions
- Screw-in baskets/paddles allow method changes in seconds
- Individual clutches allow each shaft to be raised, lowered and engaged independent of the drive head
- „Easy-Centre“ vessel centering system with safe locking and glass protection
- One piece PETG water bath with built- in drain tap: no leaks; easy to remove and clean
- Full printed test report on completion of the run
- Various options and accessories available
- Connectivity to Sampilio X8 (administration via “Browser-Interface”) *lockout spec*
- “easyDISS TX8cool” available for nano formulation – cooling down to 5°C possible *lockout spec*



easyDISS TX8

features and advantages (2)



Interchangeable Baskets/Paddles

- for method changes in seconds
- laser numbered and certificated (CoC) as standard



„Easy-Centre“ Vessel Location

- ensures that the vessels are perfectly centred very time
- vessel easy to remove and locate



one piece wather bath

- easy to remove and clean with:
 - built in drain drop
 - round edges
 - no leaks

easyDISS TX8

features and advantages (3)



easy administration for 21 CFR compliance:

If the easyDISS is connected to the Sampilio X8, it is possible to


- administrate your easyDISS with the Browser-Interface for the Sampilio easily via your (Internet-) Browser without any further software-installation and software-validation, to
- administrate easily methods, reports incl. electronic signature,
- administrate user permissions by individual user roles,
- audittrail review,
- data backup,
- parent-child-coupling,
- etc.





Teflon coated baskets and paddles

- for aggressive media

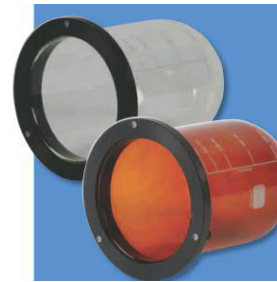

different vessels and vessel covers

- vessels with clear glass or amber
- EMC Ultra Precision Vessels
- low evaporation covers



rotating cylinders

- for transdermal patches

easyDISS TX8

options and accessories (2)




paddle over disk

- different variants available



Automatic Tablet Drop

- for simultaneous drop of dosage forms at all vessels

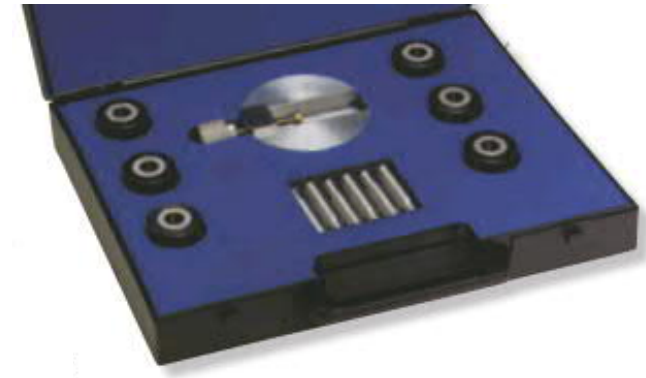

Intrinsic Kit

- easy to remove and clean with:
 - built in drain drop
 - round edges
 - no leaks



sampling accessories

- for all kinds of automatic or manual sampling



Small Volume Conversion Kit

- for special dosage forms with lower concentrations



Immersion Cell

- for testing semisolids

easyDISS TX8

options and accessories (4)



calibration equipment

- for calibration issues



accessories

- we have a big range of dissolution accessories for all brands and we will support you finding the correct ones

BI-DPX Administration / Data Management BI-SX8



Media Preparation



Dissolution Test



*Dissolution
Autosampling*

RIGGTEK

main products – customers worldwide





Thank you for your attention!
Do you have any questions?

- Visit our webpage www.riggtek.com or
- give us a call **Tel.: +49 89 2302469-0** or
- email us to support@riggtek.com

We are available!
your RIGGTEK-Team

